

Data Sheet

Model: TSH15/25/610

1500 °C Maximum Temperature,
Single Zone Horizontal Tube Furnace




INTRODUCTION

The TSH15 furnace is a bench mounted tube furnace ideal for most general laboratory thermal processing applications.

The TSH15/25 Single Zone Tube furnace can accept a work tube with internal tube diameter of 25 mm. TSH15/25 range offers a choice of three heated zone lengths of 180 mm, 450 mm & 610 mm.

TSH15/25/610 comes with a heated zone length of 610 mm.

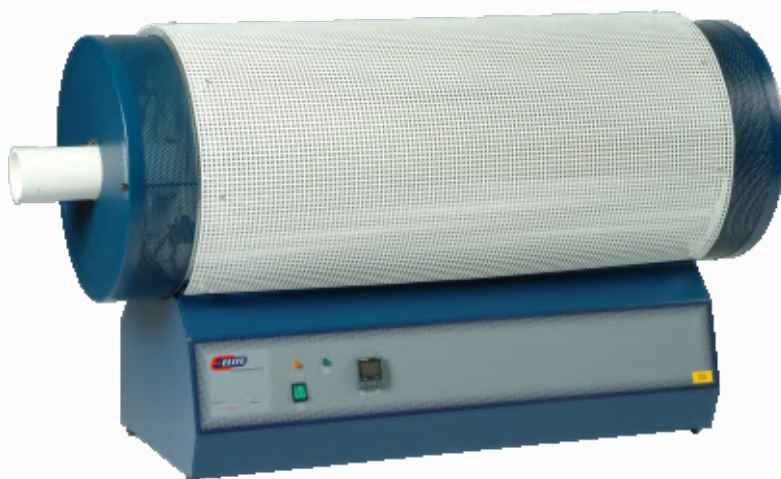

TSH15/25/610

SPECIFICATIONS

Maximum Temperature: 1500 °C

Maximum Continuous Temperature: 1500 °C

| These models are heated by axially mounted silicon carbide elements around the work tube for fast heat up and excellent temperature uniformity.

| This furnace design requires the use of a separate work tube of a grade suitable for the maximum temperature rating of 1500 °C.

| A work tube is not supplied as an integral part of the furnace and therefore needs to be ordered with the furnace as it is an essential accessory.

| Temperature Sensor: 'R' type Thermocouple.

| Energy efficient, high quality, low thermal mass insulation.

| High end Microprocessor PID controller.

| Horizontal models are supplied as standard with controls in the base.

| **External Dimensions (mm):** H x W x D **Net Wt.**
(Nominal) TSH15/25/610 715 x 1150 x 460 75 kg

| Supply: 230V– 1 Phase – 50/60Hz

| Energy rating: 6.0 kW (Average)
10.1 kW (Peak)

OPTIONAL FEATURES

| Work tubes of various materials, lengths and diameters for use in the furnace.

| The work tubes are available for containment of atmosphere or protection against process contaminants.

| A variety of triple flange gas tight end seals for work tubes to allow processing under vacuum or gas atmospheres.

| Over Temperature protection controller.

| Multi segment, multi program storage controllers.

| Thermal Radiation Plugs.

Vertical versions are available which have a separate console

Elite Thermal Systems Limited

Elite Court, 6 Stuart Road, Market Harborough, Leicestershire LE16 9PQ, UK

Tel: +44 (0)1858 469834 | E-mail: contact@elitefurnaces.com | Website: www.elitefurnaces.com